


## Multicentric ECG studies can now be performed without telephone lines.

CardiaBase and Mortara have validated a fully GSM solution for ECG centralization in multicentric, multi-countries Phase II-III studies.

This solution is based on ELI 150 GSM recorders that perform 12-lead ECG recording and automatically transmit ECG files (in native XML format) to a study server at CardiaBase, using the local GSM network. This solution requires just a basic local GSM network and works in a large number of countries.

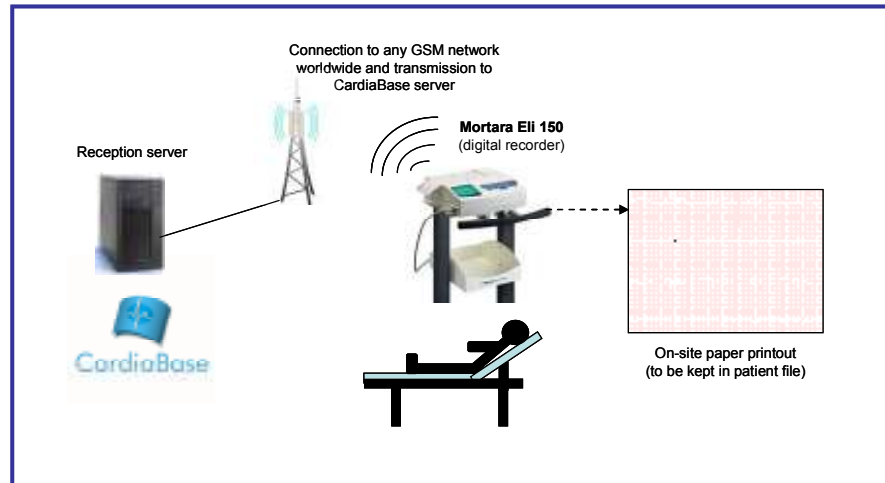

This system allows a very efficient management of sites by CardiaBase during the study, as machines:

- Are 100% programmed before shipment with study, site and data capture references (avoiding time-consuming site set up)
- Automatically roam for GSM network in any country
- Are able to transmit ECG or a full device auto-diagnosis to CardiaBase in case of problem (i.e. CardiaBase has instantaneous remote diagnosis of the machine in case of need).